


# Brochure

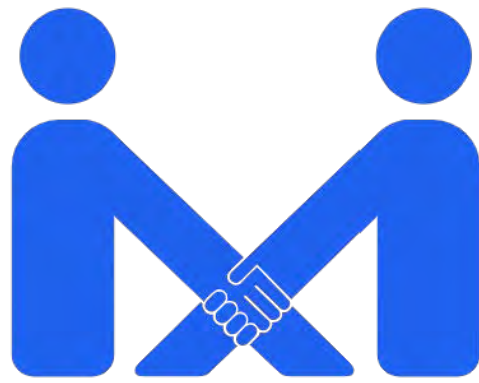

**ixicare**

# What is ixicare?


ixicare is a remote monitoring distress solution using wireless connectivity and mobile technology with indoor and outdoor location tracking, geofencing...



## Our Alert Solutions

### Manual Alerts

The manual alerts allow users to reach out when they are in need. Users are able to trigger an alert manually through the solution's wearable or immobile devices to request assistance in emergencies.

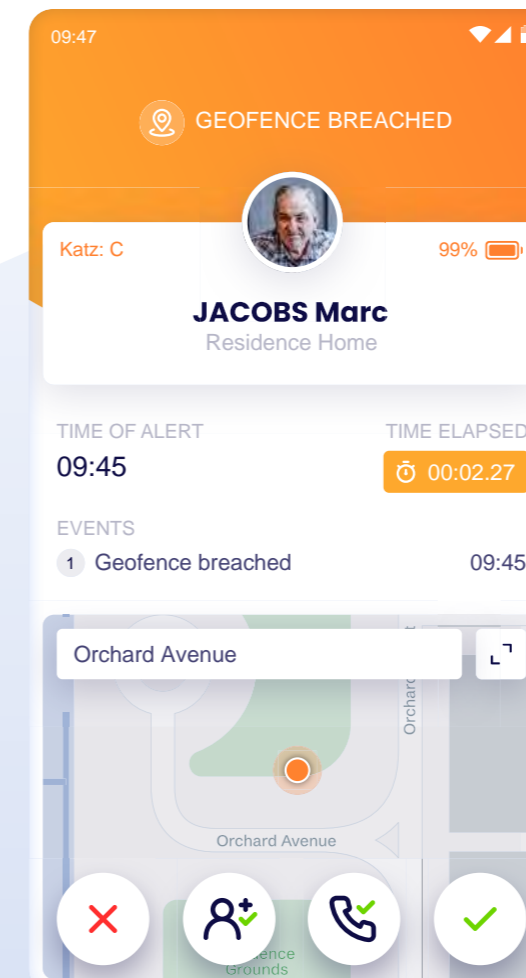

### Automatic Alerts


ixicare automatically triggers alerts through indoor and outdoor geofencing, as well as No-Motion, Man-Down, and Fall Detection. If a user cannot activate an alert manually, the system will do it automatically. Responders instantly receive the user's location on the ixicare mobile app, with optional two-way SIM voice communication for direct support.

# Our Alert Software

## ixicare Mobile Application

The ixicare mobile application can be easily downloaded from Google Play Store or from Apple's App Store.

Care providers are able to receive alerts through the mobile app. Here, they can also view the user's location.

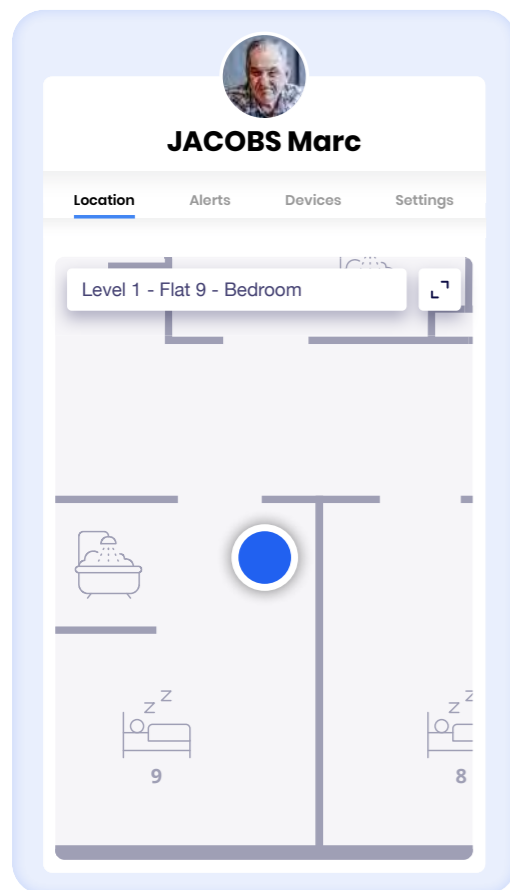


## ixicare Web Application

The ixicare web application can be accessed on any desktop computer, giving alert responders and management teams a clear overview of alerts, responses, and user activity.

Settings can be customised to define safe zones and restricted areas, while detailed reports show how much time a user spends in specific places such as the bathroom, bedroom, or kitchen. This makes it easier to spot unusual patterns that may indicate health concerns, including digestive issues, sleep disturbances, or reduced social interaction.

Beyond these features, the online web dashboard offers additional tools and insights to further support monitoring and care management.

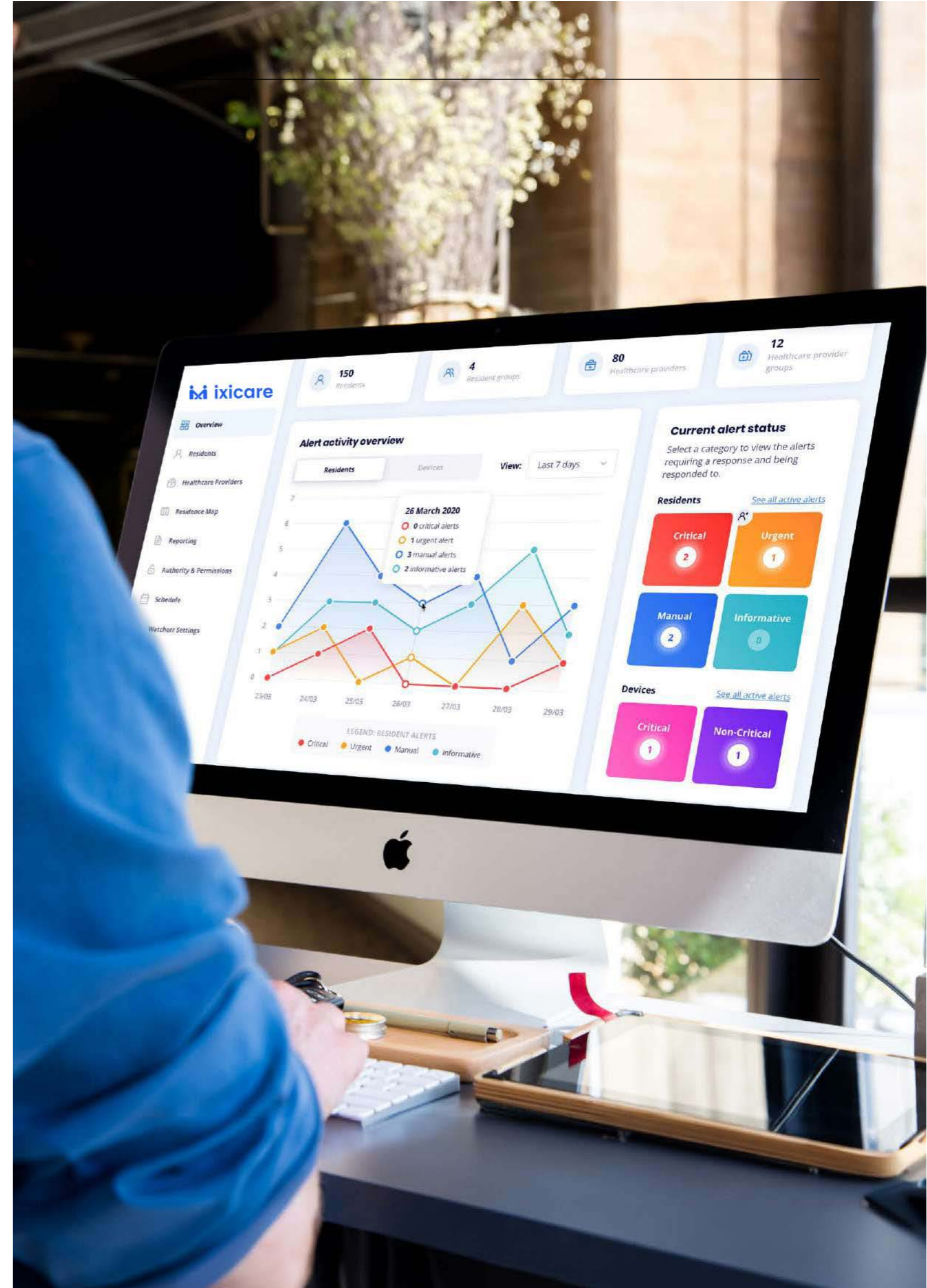

# The ixicare Devices

ixi Lite



ixi COTG



ixi Pro



ixi Touch



ixi Gateway



ixi Button



ixi Beacon

## The ixi Pro

Revolutionise your system with the ixi Pro

The ixi Pro is IP68 water-resistant, can be charged while wearing, and lets users trigger alerts manually or automatically.



### ixi Pro Features



Clock



Works indoors and outdoors



Indoor location tracking with beacons



Manual alerts via a button



Automatic geofencing alerts inside and outside via polygon



2-way SIM voice communication during an alert



Waterproof IP68



Charge while wearing



Multiple strap options



Long battery life



Automatic system updates



Fall detection, Man-down, No-motion

## The ixi Lite

\*Requires an ixi Gateway to function



### ixi Lite Features



Manual alerts with a touch of a button



Indoor location tracking with beacons



Automatic updates without the need of intervention



Battery life of more than half a year



Can be worn on the wrist or as a pendant



Waterproof

## The Ixi Button

\*Requires an ixi Gateway to function



### ixi Emergency and Alarm Button Features



No wiring or cables



Waterproof



Battery life of more than a year



Automatic system updates

## Gateways and Beacons

The ixicare solution uses Gateways and Beacons to communicate between devices.



### ixi Gateway Features



Indoor hub for all devices



Powered by PoE or a power cable



Provides a secure alert & monitoring system through ethernet cable



Reduces battery consumption of the ixi Pro



Over the air System Updates

### ixi Beacon Features




Optional: used for indoor location tracking



Wireless installation



Battery life of 9 months

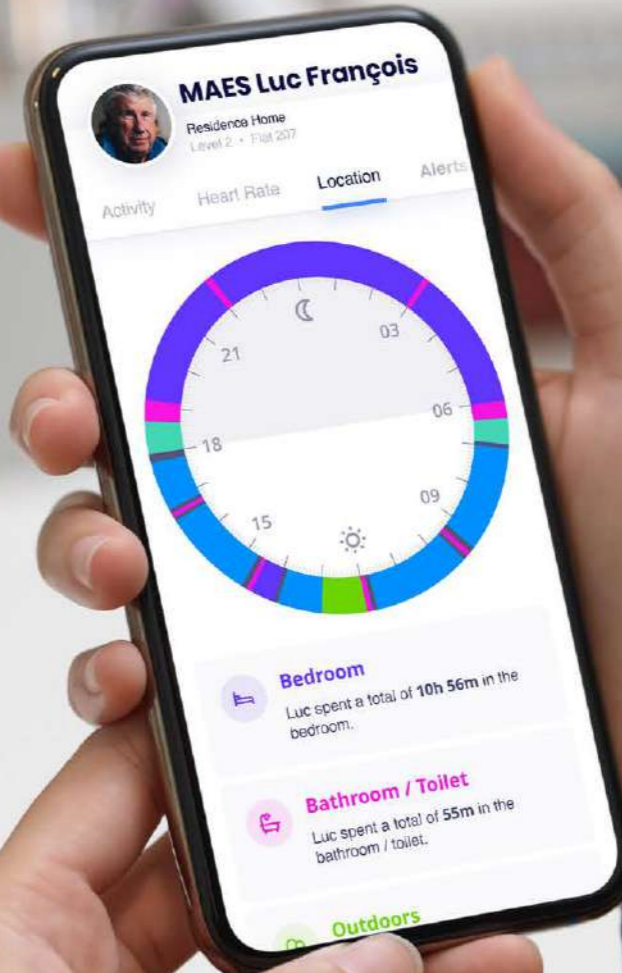


# How the system works

## Manual Alert

Our manual alert system enables users to reach out when they are in need with the press of a button. Alert responders receive the alert on their smartphone via the ixicare mobile application.

The alert contains information such as the user's location when an alert is created.

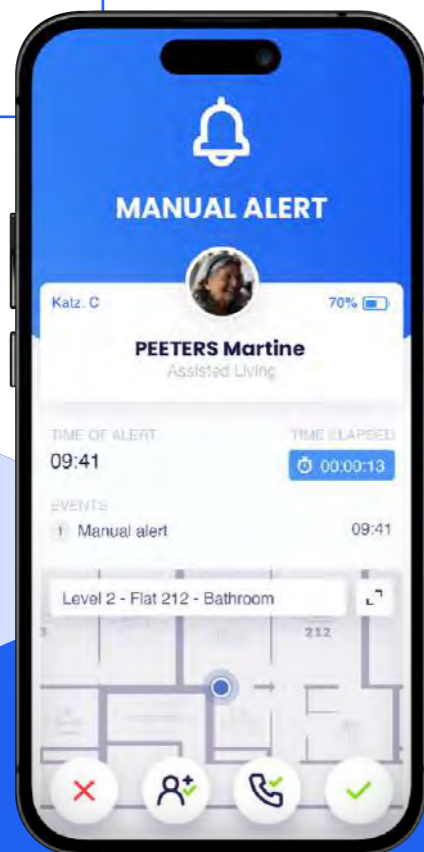

## Automatic Fall detection alert

When a user falls, an impact will be detected by the ixicare Pro, launching an automatic fall detection alert to alert responders on their ixicare mobile application.




Geofenced zones can be pre-selected on the Web App.

## Automatic GeofencingAlert

### Indoor Geofencing

Automatic indoor alerts are triggered when the user exits a predefined zone of the building. The system sends out alerts to alert responders on their ixicare mobile application. The areas can be pre-set via the web application.

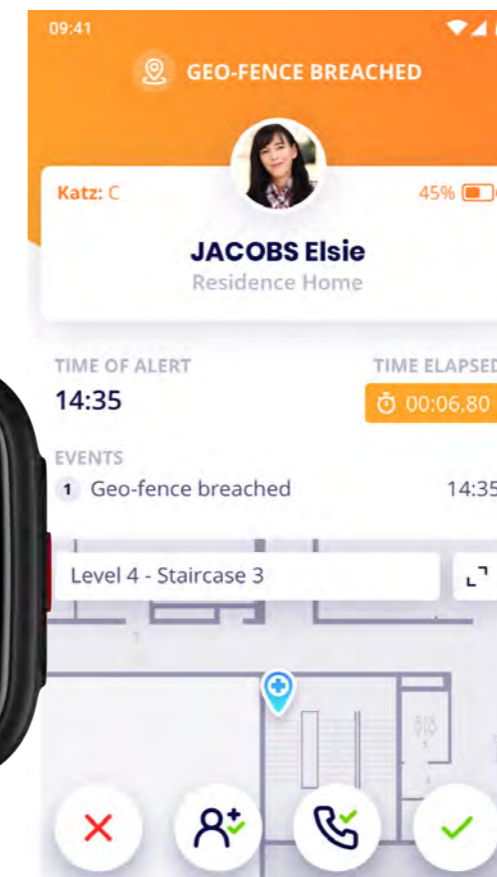

### Outdoor Geofencing

If a user exits a predefined area drawn by a polygon, an alert notification will be sent to the alert responder informing them that the user has breached the geofencing area. The notification will contain user details and a map showing where the user is. If the user exits the premises, the Pro will connect via the mobile network and activate the GPS. This way, the user can always be located and called.

# Why ixicare?

## Easy to install



Minimal Cabling Required



Wireless Buttons

## Flexible Solutions



Pro



Lite




Button

## Easy to use



Charge the Pro on the Go

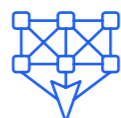

Evolves with your needs



Touch of a button



Automated alerts

## Waterproof IP68



Shower



Bath



Dishes



Swim

# ixicare Device Comparison



Feature	Pro	Lite	Touch	Button
Indoor manual alerts	✓	✓	✓	✓
Online monitoring	✓	✓	✓	✓
Waterproof	✓	✓	✓	✓
Over-the-air updates	✓	✓	✓	✓
Detect activity patterns	✓	✓	✓	✗
Indoor location tracking	✓	✓	✓	✗
Indoor geofencing alerts	✓	✓	✓	✗
Outdoor geofencing alerts	✓	✗	✗	✗
Outdoor location tracking	✓	✗	✗	✗
Outdoor manual alerts	✓	✗	✗	✗
Fall detection, Man-down, No-motion	✓	✗	✗	✗
2-way SIM voice communication	✓	✗	✗	✗
Wrist rotation activates screen	✓	✗	✗	✗
Charge while wearing	✓	✗	✗	✗
Clock	✓	✗	✗	✗
Touch screen	✓	✗	✓	✗
Battery Life*	+/- 7 days	+/- 6 months	+/- 6 months	+/- 1 year


\* Depends on activated features

# Easy plug and play solution

Our cable-free, economical solution is perfect for both new buildings and renovations, can be integrated in your system, and can be offered via private label. We provide guidance and training so your staff can easily manage it day-to-day.




ixicare can be used as a standalone system or be integrated in your systems

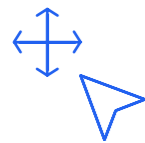

Add, remove and move devices




Monitor device performance and battery levels

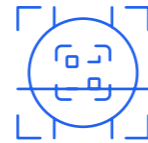

Scan lost devices to identify to whom they belong.



# Infrastructure

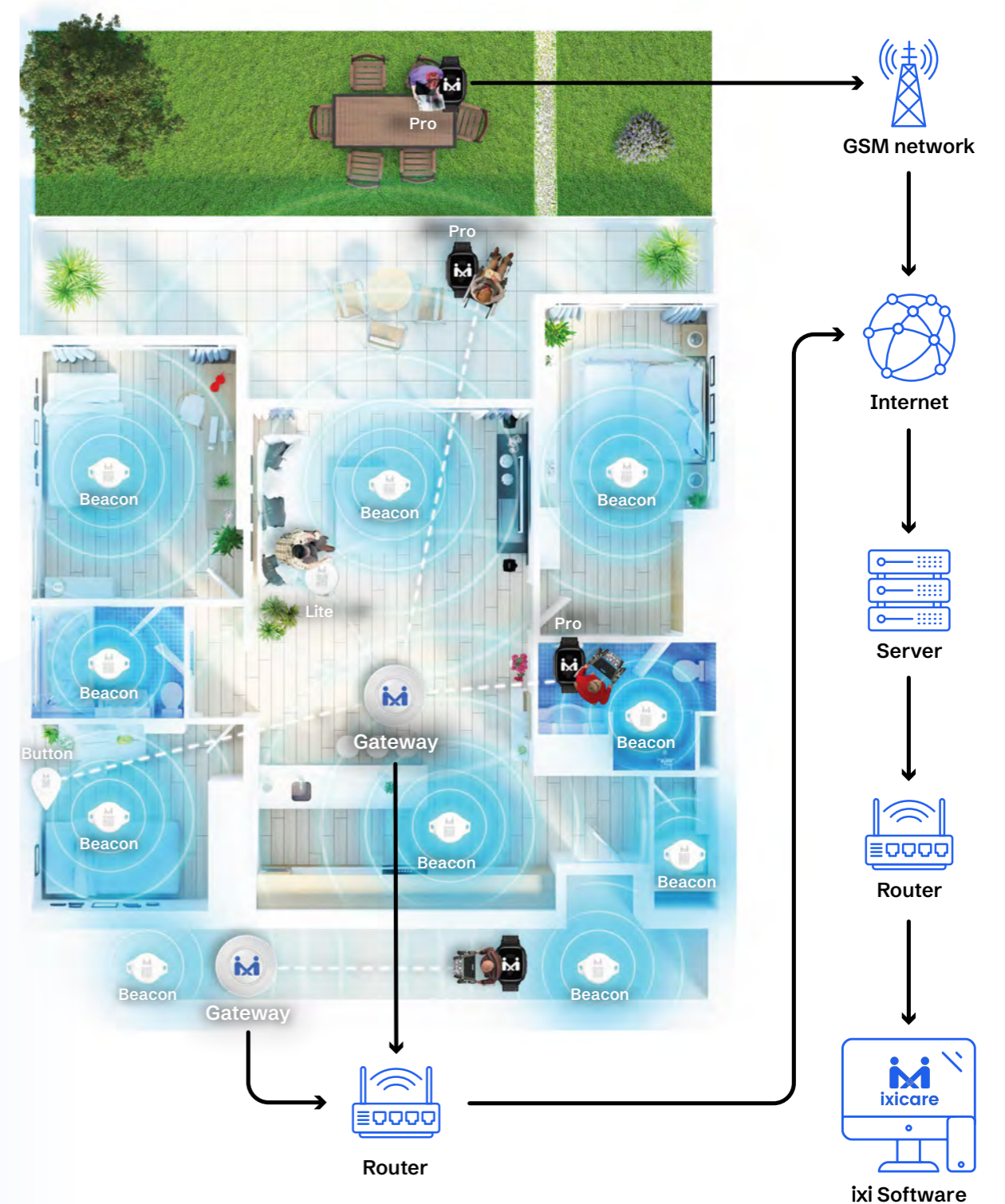



Copyright © 2025 ixicare  
BV All rights reserved.

No part of this publication may be reproduced, distributed, or transmitted in any form or by any means, including photocopying, recording, or other electronic or mechanical methods, without the prior written permission of ixicare.

For more information, contact us on  
[info@ixicare.com](mailto:info@ixicare.com)

# Privacy & Security



ixicare is guided and supported by a specialised GDPR agency



We strive to be secure, straightforward and up-to-date in our data practices



We will not access personal data without your explicit request to do so



Secure cloud data storage on European servers



## Contact

[support@ixicare.com](mailto:support@ixicare.com)

Warande 13B

9270 Laarne, Belgium

[www.ixicare.com](http://www.ixicare.com)